

DrägerService®

Dräger

NARKOMED 6000 TECHNICAL SERVICE BULLETIN # 6001

Product: Narkomed 6000

Date: October 30, 2000

Re: Power cord routing

Reference Doc: ERN # 00-0783

Reason/Rational: Document proper power cord routing to prevent unintended activation or deactivation of circuit breakers.

Solution/Procedure: Route the power cord's strain relief with a loop in a clockwise direction. The retaining clip shall to be located $16.5 \pm .250$ inches, (5.25 inch diameter loop approx.) from the cord end and secured to the Vitalbus HUB door using the mounting hole directly above the power cord socket. See attached illustration, Figure - 1. Refer to Narkomed 6000 Technical Service manual Replacement Procedures, Section - 4 for complete instructions.

Devices affected: All Narkomed 6000 shipped before November 2nd 2000.

Timing/When: Next Service Call, Periodic Manufacturer's Certification, or Vigilance Audit®.

Additional Info: n/a

Cost: n/a

Ordering Info/Parts: n/a

Distribution: All Authorized Service Organizations, DrägerService

If you have any questions, please contact Technical Support by phone at 1-800-543-5047 or by fax at 1-215-721-5789.

DrägerService®

Dräger

NARKOMED 6000
TECHNICAL SERVICE BULLETIN # 6001

Page - 2

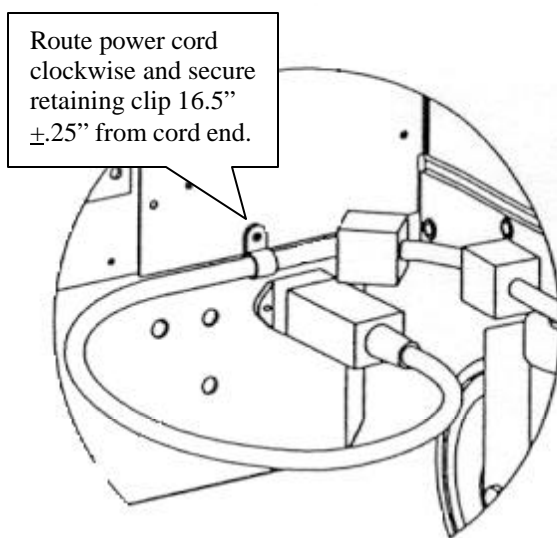


Figure - 1 Power Cord Routing Detail